


# Portable Ventilator

**SIARO 300** 





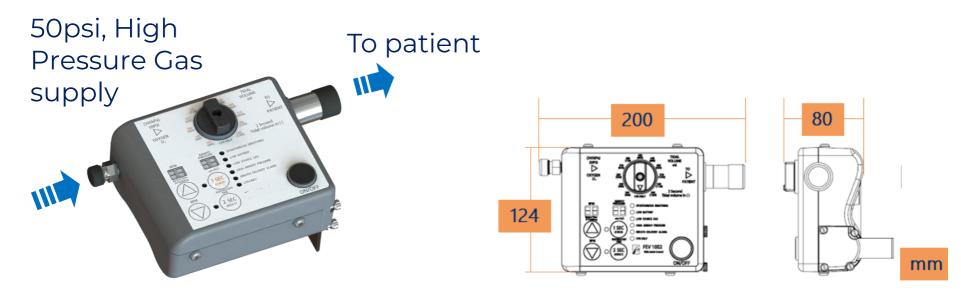


- 1. 100% Oxygen Gas powered, Electronically controlled ventilator.
- 2. <u>Durable and Lightweight with 2 D size batteries. < than 1Kg.</u>
- 3. Patient type: Bodyweight greater than 20 kg.



- 1. Trigger sensitivity at less than 2cmH<sub>2</sub>O to assist spontaneous breathing.
- 2. Independent setting of inspiratory time, tidal volume and BPM control.
- 3. Pre-set visual and audio alarm.
- 4. "CPR only" Special mode design.

## Replacing of Batteries

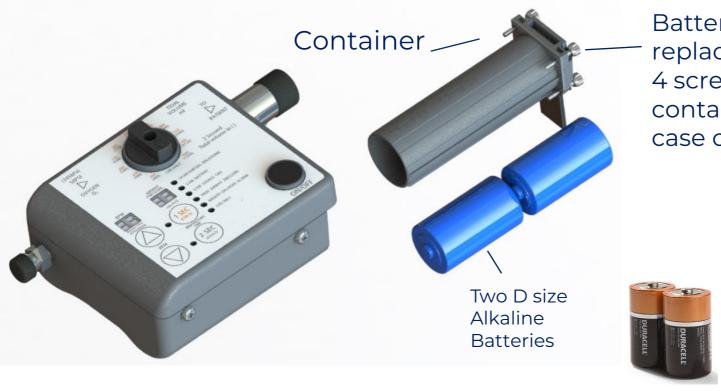

Batteries can be replaced by removing 4 screws and pull the container out of the case of the ventilator.




## **Specialty Product**

Streamlined design on the back of the case, easy to hold



IPX4 waterproof, Protecting from splashing water

## Model Conmed Photo

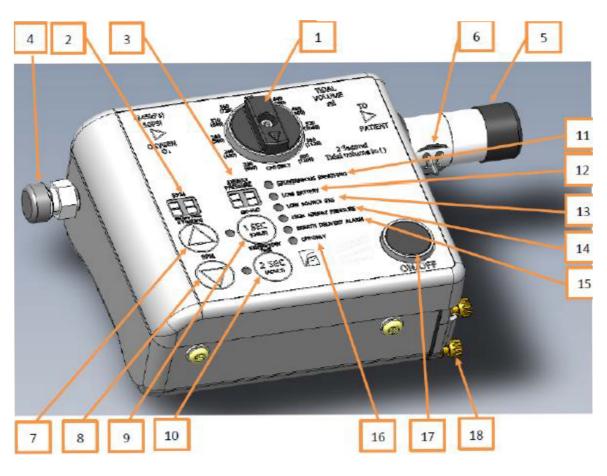

### 3) Components:

- 1 Tidal volume control
- 2- BPM Display
- 3- Airway pressure display
- 4- Oxygen inlet
- 5 -Patient circuit connection
- 6- Anti-asphyxia valve
- 7- BPM control (increase)
- 8- BPM control (decrease)
- 9- Inspiratory time (1 second) control

- 10- Inspiratory time (2 seconds) control
- 11-Spontaneous breathing indicator
- 12- Low battery alarm
- 13- Low source gas alarm
- 14- High airway pressure alarm
- 15- Breath delivery alarm
- 16- CPR only mode
- 17- Power ON/OFF switch
- 18- Battery replacement lock screw

accessory:





CONTACTS,

S.S Padana Superiore, 317 20055

Vimodrone (MI) ITALY

http://www.ventexmedical.com

info@ventexmedical.com